

AVImark ReadMe 2016.2.2

Table of Contents

Please use the bookmark menu on the left to navigate between each project and feature.

- 1. Patient Specific Reminder Enhancements**
- 2. Enhancement to E-mail Directly from Invoice Window**
- 3. 2016.2.2 Fixes List**

Patient Specific Reminder Enhancement

Enhancements were made to patient specific reminders in the AVImark software so they are now able to interact as a member of a reminder series. The work flow of patient specific reminders was modified to now compare the medical history date of existing series members to the modified date of the patient specific reminder and either update, replace, or remove the reminder based on the most recent modified date.

In addition to the modifications of the workflow, there are several other areas of reminders that were enhanced.

Must choose a Valid Code – when entering a code in the patient specific reminder field it now requires a code that exists in either the Treatment, Inventory, Diagnosis, or Problem list. It will no longer allow a code to be entered that does not exist within the program.

AVI
mark

Patient Reminder

Reminder

REMINDER

Code [] Description []

Due []

Choose

- Treatment F7
- Inventory Shift+F7
- Diagnosis Ctrl+F7
- Problems Shift+Ctrl+F7


Controls Advanced Reminded By

Remind every 1 Year(s) Remind if never treated

Start reminding At Age [] (none) []

Stop reminding At Age [] (none) []


If a code that does not exist in any of the lists above is entered into the code field the user will receive an alert that an invalid code has been entered and to please enter a valid code.

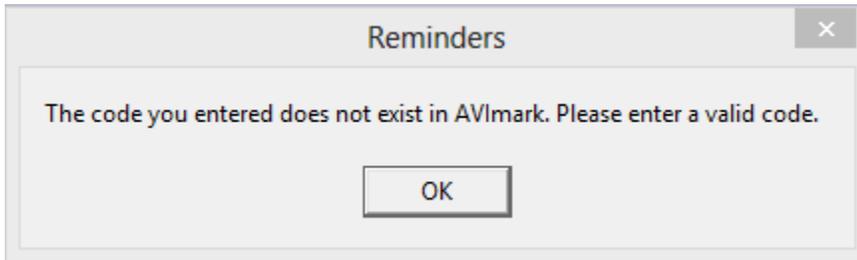

Locked Fields in Patient Specific Reminder Window - When entering a new patient specific reminder some of the fields within the reminder template window will be grayed out depending on if the reminder code selected has a pre-existing template setup on it or not.

Reminder Code with a Template:

If the selected code has a pre-existing reminder template the only fields available to modify are the Due date on the controls tab or Priority, Suspend, and Resume After fields on the advanced tab.

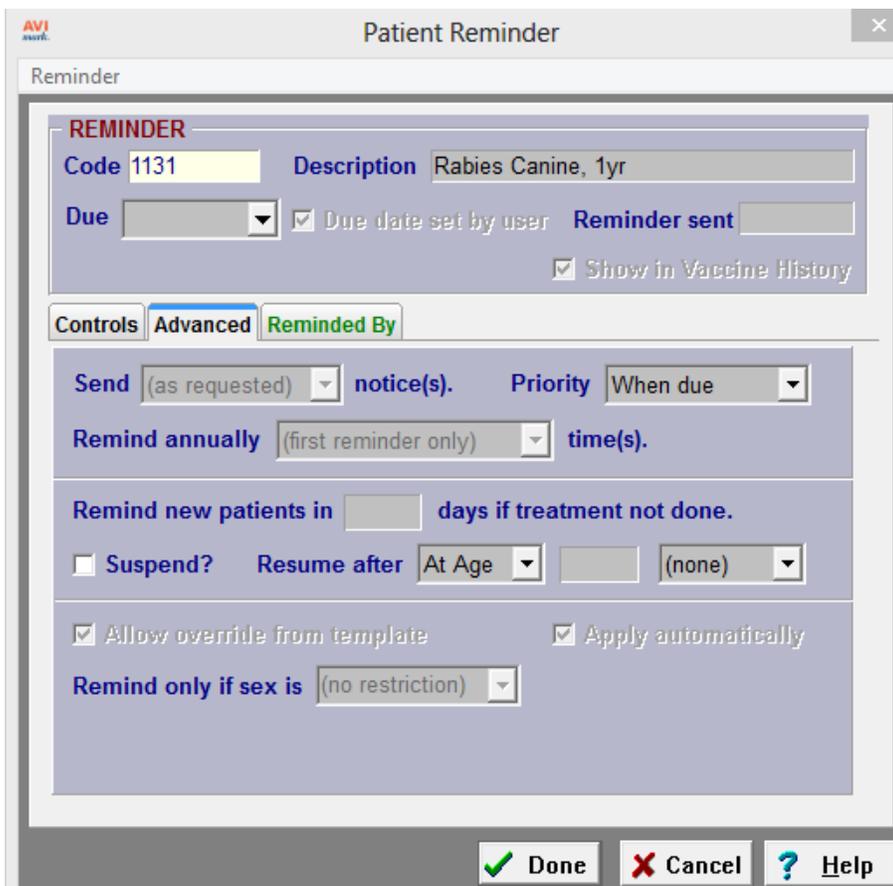

Reminder Code without a Template:


If the selected code does not have a pre-existing reminder template associated to it the user can modify the Due date, Remind every, and Patient to be treated fields on the controls tab or Priority, Suspend, and Resume After fields on the advanced tab.

The screenshot shows the 'Patient Reminder' dialog box with the 'Controls' tab selected. The 'REMINDER' section includes fields for 'Code' (DENTEX), 'Description' (Dental Exam), 'Due' (dropdown), and 'Reminder sent' (text field). A checkbox for 'Show in Vaccine History' is checked. Below this are three tabs: 'Controls', 'Advanced', and 'Reminded By'. The 'Controls' tab contains: 'Remind every' (1) 'Year(s)' (dropdown), 'Remind if never treated' (checkbox), 'Start reminding' (At Age dropdown, text field, (none) dropdown), 'Stop reminding' (At Age dropdown, text field, (none) dropdown), 'Patient to be treated' (text field, time(s)), and 'Series name' (text field with up/down arrows).

The screenshot shows the 'Patient Reminder' dialog box with the 'Advanced' tab selected. The 'REMINDER' section is identical to the Controls tab. The 'Advanced' tab contains: 'Send' (as requested dropdown), 'notice(s)' (text field), 'Priority' (When due dropdown), 'Remind annually' (first reminder only dropdown), 'time(s)' (text field), 'Remind new patients in' (text field) 'days if treatment not done.', 'Suspend?' (checkbox), 'Resume after' (At Age dropdown, text field, (none) dropdown), 'Allow override from template' (checkbox), 'Apply automatically' (checkbox), and 'Remind only if sex is' (no restriction dropdown).

Note: Once a patient specific reminder has been added to the patient area the only fields the user can edit is the Due date, Priority, Suspend, or Resume After fields.

Requires a Due Date – The Due date field must be entered before you can save a new patient specific reminder. If the due date is left blank the user will get an alert when clicking done that the date is blank and one must be selected before continuing.

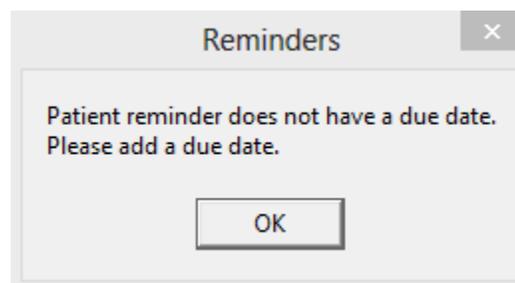

Enhancement to E-mail Directly from Invoice Window

The AVImark invoice window has been enhanced to include an e-mail button that will allow the user to e-mail the customer a copy of their invoice directly from the invoice window.

The screenshot shows the 'Invoice' window for AVImark Animal Hospital. The window title is 'Invoice' and it has a menu bar with 'Invoice', 'Tools', and 'Insurance'. The main content area is divided into several sections:

- Date:** 05-02-16
- AVImark Animal Hospital**
- Invoice no. (open)**
- Sold To:** John Doe, 1234 Smith Street, Columbia, SC 29201
- SSN:** (blank), **D/L:** (blank)
- Patients:** A table with one entry: 'Test' with an amount of 0.00.
- Table of Services:**

Date	Patient	Code	Description	Qty	Amount
05-02-16	Test	7550	Extract Tooth	0	0.00
- Summary:** Balance: 0.00, Charges: 0.00, Tax: 0.00, Payments: 0.00, New balance: 0.00
- Send To:** OneNote 2013 (dropdown), **Picture:** (dropdown), **Copies:** 1 (input), **Preview:** (checkbox), **Include notes:** (checkbox)
- Buttons:** Print, Email, Close, Quit, Help

Post copy to Medical History – within the e-mail window there is also an option to have an entry posted to the patient’s medical history for record tracking purposes.

The screenshot shows the 'Send e-mail' window. The window title is 'Send e-mail'. The main content area contains the following fields and options:

- Send To:** johndoe@fakeaddress.com
- CC:** (empty)
- Subject:** (empty)
- Attachment:** C:\Users\JDoe\AppData\Local\Temp\AVIDAD48DE9.PDF
- Reply To:** owa.avimark.net (dropdown)
- Post copy to medical history:** (checked checkbox)
- Patient:** Test (dropdown)
- Send in plain text:** (selected radio button), **Send as HTML content:** (radio button)
- Buttons:** Print, Ricoh Aficio MP 4500 PCL6 (dropdown), Send, Quit

2016.2.2 Fixes List

7718 – Fix Memory Leak Updating Inventory Sales Statistics

Updating inventory sales statistics would display a memory leak error in large databases.

7703 – Fix Ability to Add Unprocessed Payments Without it Adding a Partial Payment

If an unprocessed payment was added into accounting on a computer that did not process credit cards it would go in as a partial payment.

7720 – Add Right Click “Re-Enter” Option to Replace “Refill” Option on Services in Medical History

The “re-enter” option was added to services in medical history to replace the “refill” option so clinics still have the ability for the service to go back into history and carry over previous label notes.

7789 – Fix Reminder Merge Word Displaying “No Reminders Due”

The <reminders> merge word would display “No reminders due” on documents when the patient actually had current reminders.

7796 – Optimize Period Totals

Changes in the way the Period Totals report calculated new clients caused increased slowness running the report.

2016.2.1 Fixes List

5011 – Include Inventory Used When Linking a Canned Estimate to an Estimate

When a canned estimate containing a service that included inventory used was added to a patient specific estimate, the inventory used was not transferring over.

7495 – Access Violations Printing Treatment Labels

Printing a treatment label either directly from the treatment list or from the entry in medical history would result in an access violation.

7496 - Fix Treatments So Saved Labels Carry Over Into Medical History

Labels attached to a treatment would not carry over into medical once that treatment was added into medical history.

6988 - Axis-Q Results Should Still Auto-Import After HS Sync Services Are Restarted

If the Axis-Q HS Sync services were restarted the Axis-Q results would not auto-import until the Axis-Q sync service was stopped and restarted on the server as well.

2208 - Scriptel Signature Too Small to View on Documents

When viewing document signatures from an attachment area they appeared very small or non-existent.

2016.2.0 Fixes List

6924 – Word Document Locks When Accessed By Multiple Computers across a Network

If more than one person would try to print a MS Word document at the same time across the network the document would lock and cause errors.

615 – Scriptel Device Not Prompting For Signature if Document is Previewed

The Scriptel device would not prompt for a signature on a MS Word document unless the document was printed and not if it was just previewed.